
Electronic Supplementary Material (ESI) for Lab on a Chip. This journal is © The Royal Society of Chemistry 2020

- 1 Functional angiogenesis requires microenvironmental cues balancing endothelial cell migration and proliferation
- 2 William Y. Wang¹, Daphne Lin¹, Evan H. Jarman¹, William J. Polacheck², and Brendon M. Baker^{1*}
- 3 ¹Department of Biomedical Engineering, University of Michigan, Ann Arbor, MI 48109, USA
- 4 ²Joint Department of Biomedical Engineering, University of North Carolina at Chapel Hill and North Carolina State
- 5 University, Chapel Hill, NC 27514, USA
- 6 *Correspondence to: bambren@umich.edu
- 7

8 SUPPLEMENTARY MATERIAL

9 The Supplementary Material includes 6 Supplementary Figures, 1 Supplementary Table, and 5 Supplementary Movies.

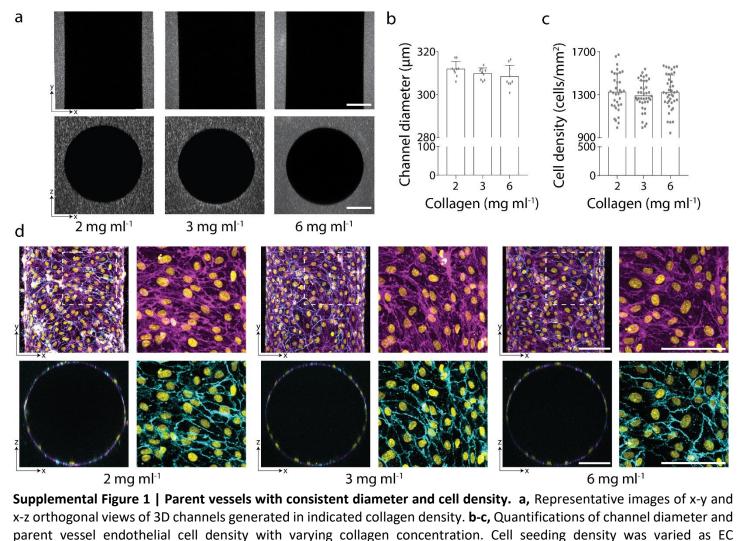


- 14 attachment efficiency was dependent on collagen density (1.5 M/ml for 2 mg ml⁻¹, 2 M/ml for 3 mg ml⁻¹ and 5 M/ml for 6
- 15 mg ml⁻¹). All data presented as mean \pm s.d.; * indicates a statically significant comparison with P<0.05 (one-way analysis 16 of variance). n≥8 channels (**b**) and n≥34 vessel segments (each 400 µm length) (**c**) per condition. **d**, Representative images
- 17 of x-y (max intensity projection) and x-z (single slice) orthogonal views of parent vessels generated with varying collagen
- 18 concentration after 24-hour culture. Insets indicated with dashed white lines. Nucleus (yellow), F-actin (magenta), and VE-
- 19 cadherin (cyan).


 $\begin{array}{c} 10 \\ 11 \end{array}$

12

13

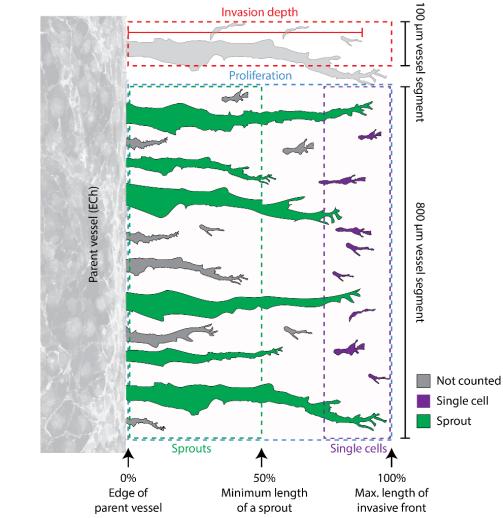


20 $\overline{2}1$ Supplemental Figure 2 | Schematic of quantification metrics. To quantify invasion depth (red), images were segmented 22 into 100 µm segments along the parent vessel and measurements were taken from the ECh edge to the tip of the 23 endothelial cell that invaded furthest within that segment. To quantify proliferation (blue), images were segmented into 24 800 μm segments along the parent vessel, and the percentage of EdU⁺ nuclei was measured only in endothelial cells that 25 had invaded into the extracellular matrix from the ECh edge to leading invasive front. To quantify leading edge single cells 26 (purple), images were segmented into 800 µm segments along the parent vessel and isolated single cells were quantified 27 within 150 µm of the leading invasive front. To quantify multicellular sprouts (green), images were segmented into 800 28 µm segments along the parent vessel, and quantified connected endothelial sprout structures (from ECh edge to sprout 29 tip) with a minimum length of half the max invasion depth. Single cells and sprouts outside of these criteria were not 30 included in the analysis (grey).

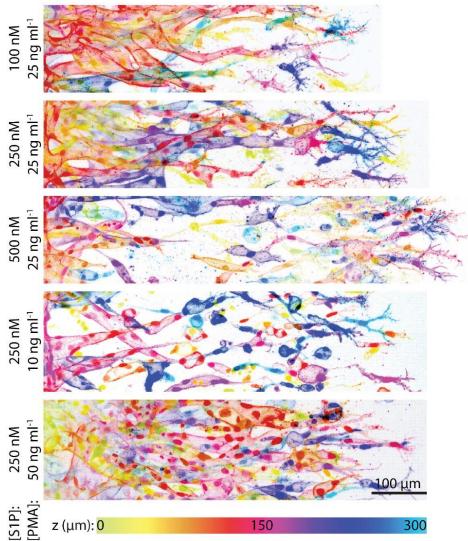


31 32

Supplemental Figure 3 | Z-depth encoded projections. 300 µm max intensity projections from Figure 2 pseudo-colored 33 with z-position. Z-stack image files were processed in FIJI using "Temporal-Color Code" function to assign each z-slice an 34 individual color. "Physics" color scale was utilized and inverted. 35

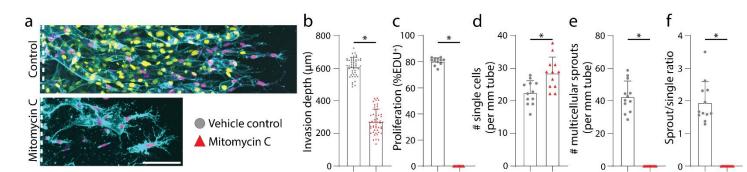


36 37

Supplemental Figure 4 | Proliferation inhibition abrogates multicellular sprouting. a, Representative images (max 38 intensity projection) of endothelial cell invasion in response to mitomycin C proliferation inhibition with S250:P50 and 39 3mg ml⁻¹ collagen. **b-f**, Quantifications of invasion depth, proliferation, and morphology of invading endothelial cells as 40 single cells or multicellular sprouts. All data presented as mean ± s.d.; * indicates a statically significant comparison with 41 P<0.05 (two-tailed Student's t-test). For invasion depth analysis (b), n=48 vessel segments (each 100 μ m length) per 42 condition. For proliferation and migration mode analysis (c-f) n=12 vessel segments (each 800 μ m length) per condition 43 from n=2 devices/condition (technical replicates) over $n \ge 4$ independent studies (biological replicates).

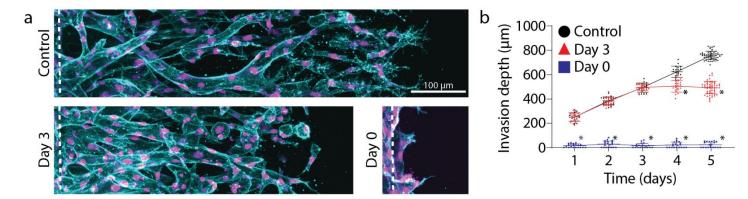

- 44 45 Supplemental Figure 5 | S1P receptor inhibition abrogates S1P-driven EC invasion. a, Representative images (max 46 intensity projection) of invading endothelial cells in response to 100 nM FTY720 treatment with S250:P25 and 3mg ml⁻¹ 47 collagen. Day 0 conditions were composed of ECs pre-treated with FTY720. F-actin (cyan), nucleus (magenta). b, 48 Endothelial cell invasion depth over time in response to 100 nM FTY720 treatment. All data presented as mean ± s.d.; * 49 indicates a statically significant comparison with P<0.05 (one-way analysis of variance). $n \ge 20$ vessel segments (each 100 50 μm length) per condition.
- 51

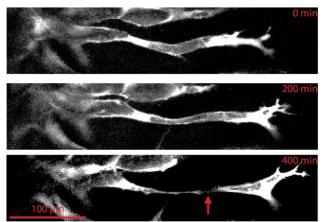

- 52 53 Supplemental Figure 6 | Tip cell breakage from stalk cell. Time-lapse series demonstrating the separation of tip and stalk
- 54 cell connectivity of a sprout cultured within 3 mg ml⁻¹ collagen with 250 nM S1P and 25 ng ml⁻¹ PMA. Time-lapse imaging
- 55 began on day 2 of sprout culture. Red arrow indicates tip and stalk cell breakage.
- 56

Figure #	[S1P] (nM)	[PMA] (ng ml ⁻¹)	[Collagen] (mg ml ⁻¹)
1с-е	0	0	3
2a	0-500	0	3
2b	0-500	25	3
2h	250	0-50	3
3a	100-500	10-50	3
4a	250	25	3
5	250-500	10-25	3
ба	250	25	2-6
6b	250	25	2-6
6i	100-250	25	2
Supp. 1d	0	0	2-6
Supp. 4	250	50	3
Supp. 5	250	25	3
Supp. 6	250	25	3

- 57 Supplemental Table 1 | Microenvironmental cues utilized in experiments by figure.
- 58
- 59

60	Supplemental Movie 1 Endothelial sprouts cultured with 100 nM S1P and 25 ng ml ⁻¹ PMA in 3 mg ml ⁻¹ collagen over 5
61	days.
62	
63	Supplemental Movie 2 Endothelial sprouts cultured with 250 nM S1P and 25 ng ml ⁻¹ PMA in 3 mg ml ⁻¹ collagen over 5
64	days.
65	
66	Supplemental Movie 3 Endothelial sprouts cultured with 500 nM S1P and 25 ng ml ⁻¹ PMA in 3 mg ml ⁻¹ collagen over 5
67	days.
68	
69	Supplemental Movie 4 Endothelial sprouts cultured with 250 nM S1P and 10 ng ml ⁻¹ PMA in 3 mg ml ⁻¹ collagen over 5
70	days.
71	
72	Supplemental Movie 5 Endothelial sprouts cultured with 250 nM S1P and 50 ng ml ⁻¹ PMA in 3 mg ml ⁻¹ collagen over 5
73	days.