

A simple, quick and non-destructive approach for sampling drugs of abuse in tablets and blotter for qualitative analysis by paper spray mass spectrometry

Mariana C. C. Diniz,^a Fabiana de Moura,^{a,c} Yuri Machado,^{*b} José Coelho Neto,^b and Evandro Piccin^a

^a Departamento de Química, Instituto de Ciências Exatas, Universidade Federal de Minas Gerais, Antônio Carlos Str., 6627, CEP 31270-901, Belo Horizonte, MG, Brazil.

^b Seção Técnica de Química e Física Legal – Divisão de Laboratório de Pesquisa, Instituto de Criminalística da Polícia Civil de Minas Gerais, Juiz de Fora Str., 400, CEP 30180-060, Belo Horizonte, MG, Brazil.

^c Departamento de Química, Centro Federal de Educação Tecnológica de Minas Gerais, Amazonas Str., 5253, CEP 30421-169, Belo Horizonte, MG, Brazil.

* yurimachad@gmail.com, +55 31 98477-9786

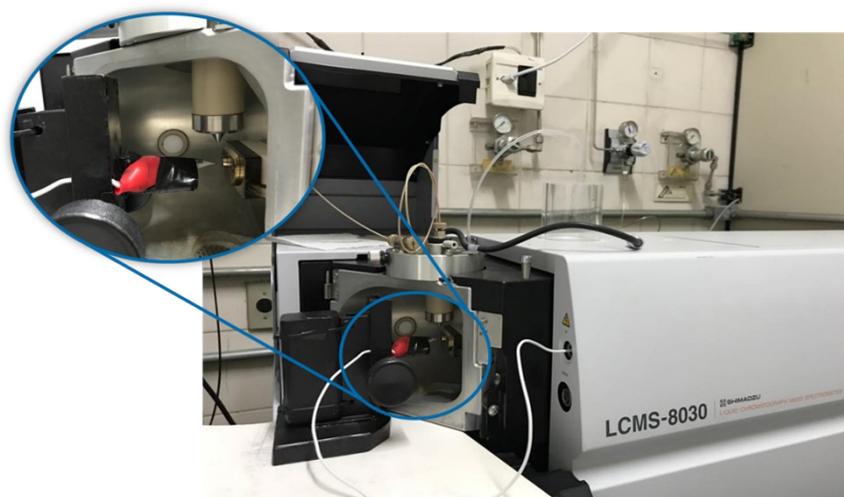


Fig. 1S Picture of the PS interface adapted on a Shimadzu mass spectrometer model LCMS 8030.

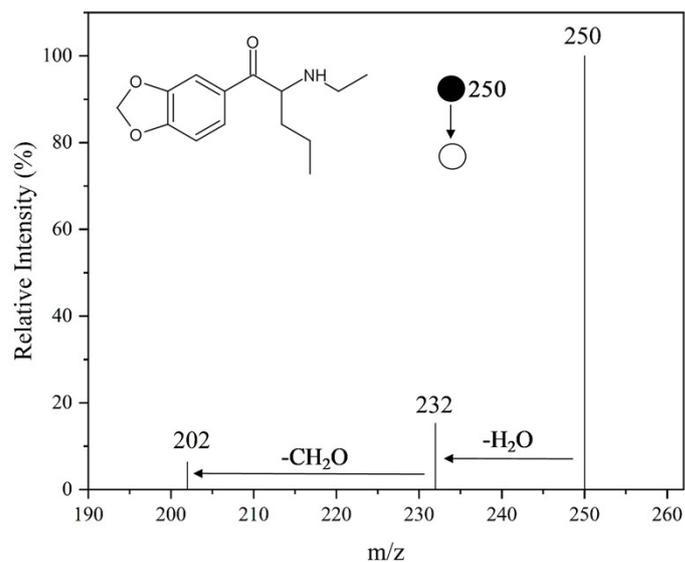


Fig. 2S PS(+)-MS/MS of *N*-ethylpentylone (NEP) recorded from a sample also containing U-47700 (a potent synthetic opioid also known as pink heroin). Conditions: support, filter paper; spray solvent, MeOH; spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision 5 V.

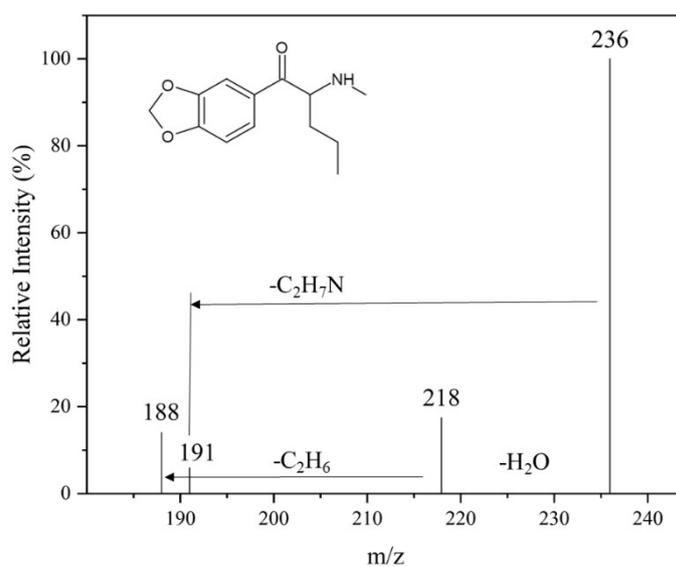


Fig. 3S PS(+)-MS/MS of pentylone recorded from a sample also containing NEP. Conditions: support, filter paper; spray solvent, MeOH; spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision 5 V.

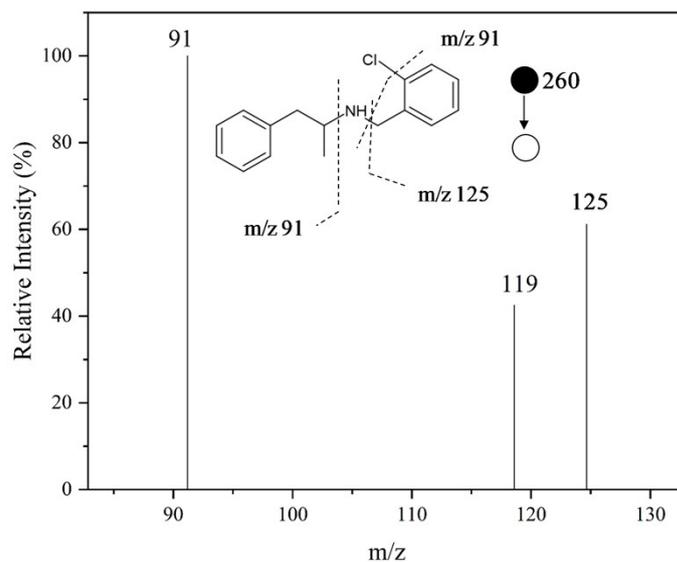


Fig. 4S PS(+)-MS/MS of clobenzorex. Conditions: support, blotter paper; spray solvent, MeOH:ACN (1:1 v/v); spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision, 5 V.

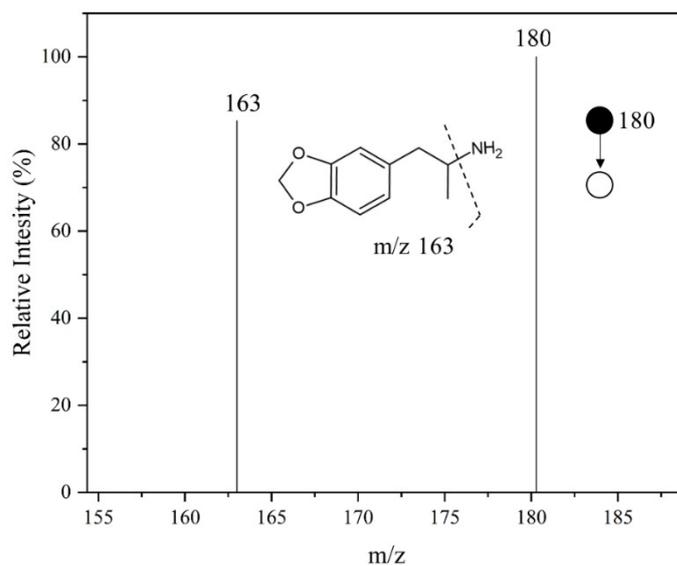


Fig. 5S PS(+)-MS/MS of MDA. Conditions: support, blotter paper; spray solvent, MeOH; spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision, 5 V.

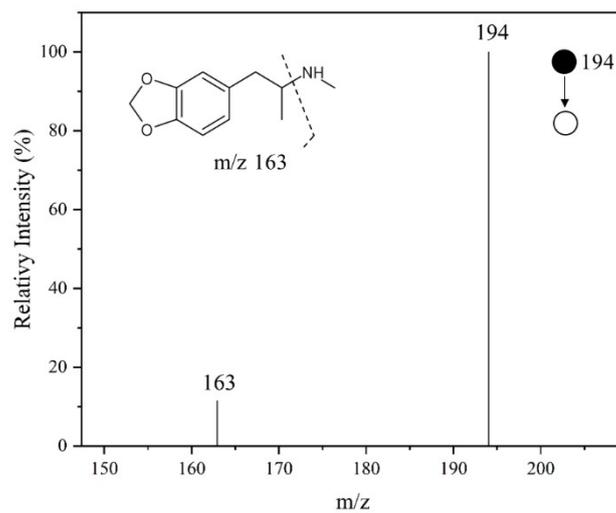


Fig. 6S PS(+)-MS/MS of MDMA. Conditions: support, blotter paper; spray solvent, MeOH; spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision, 5 V.

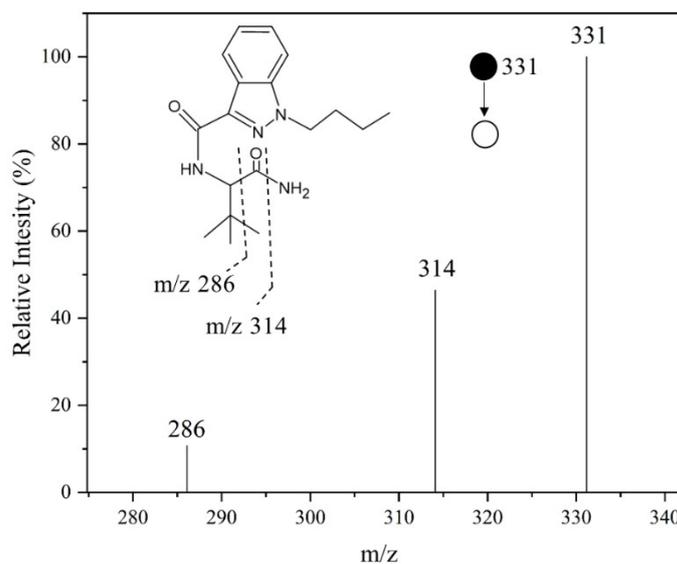


Fig. 7S PS(+)-MS/MS of ADB-BUTINACA. Conditions: support, filter paper; spray solvent, MeOH; spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision, 5 V.

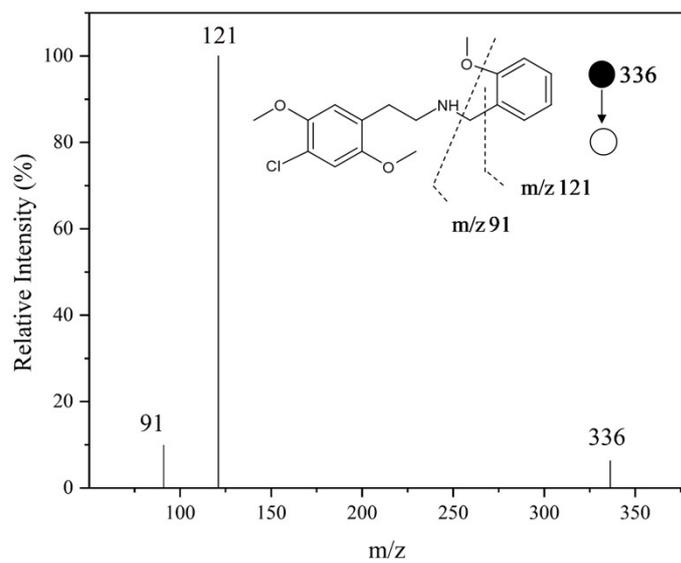


Fig. 8S PS(+)-MS/MS of 25C-NBOMe. Conditions: support, filter paper; spray solvent, MeOH; spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision, 5 V.

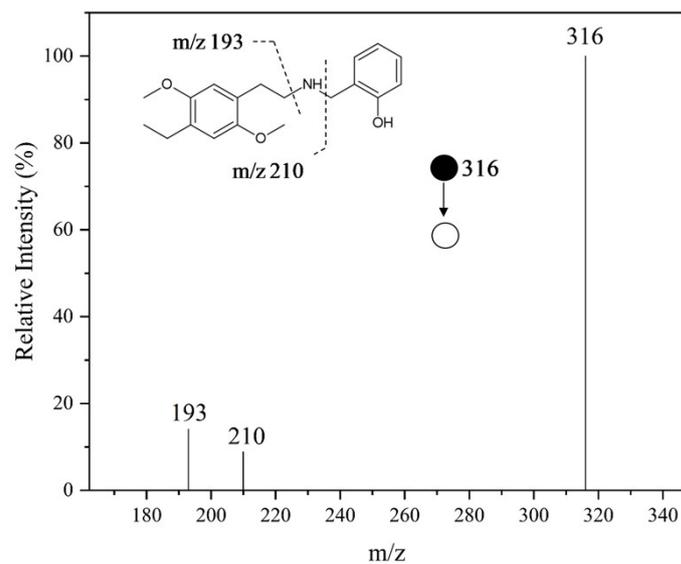


Fig. 9S PS(+)-MS/MS of 25E-NBOH. Conditions: support, blotter paper; spray solvent, MeOH:ACN (1:1 v/v); spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision, 5 V.

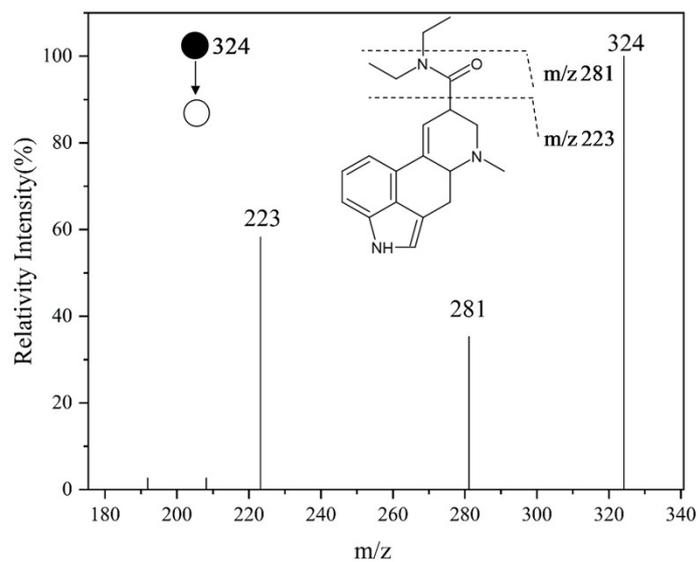


Fig. 10S PS(+)-MS/MS of 25C-NBOMe. Conditions: support, filter paper; spray solvent, MeOH; spray solvent volume, 10 μ L; spray voltage, 4.5 kV; distance between paper tip and MS inlet, 10 mm; energy collision, 20 V.

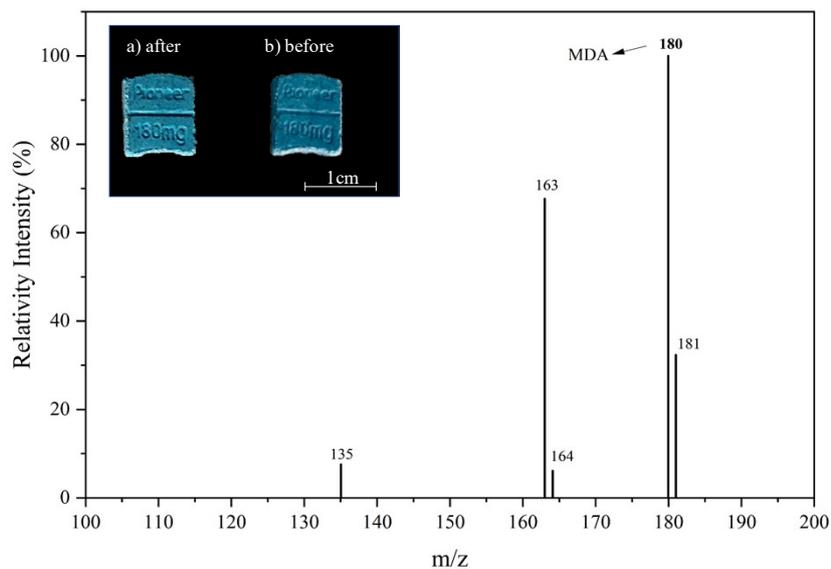


Fig. 11S PS(+)-MS spectrum of MDA after the 10 consecutive extractions with blotter paper and methanol as spray solvent. Photos of sample after 1 extraction and 10 consecutive extractions.