

## Supplementary Information

### **Upconversion Fluorescence Immunochromatographic Assay for Highly Sensitive and Rapid Detection of Ovarian Function Indicators Follicle-stimulating hormone and anti-Müllerian hormone**

Jianfeng Yao,<sup>#a</sup> Ming Li,<sup>#b</sup> Rongfu Huang,<sup>c</sup> Yanting Wang,<sup>a</sup> Liang Song,<sup>\*d</sup> Liying Chen,<sup>\*a</sup> and Zhenzhu Zheng,<sup>\*a</sup>

<sup>a</sup> Quanzhou Maternity and Child Healthcare Hospital, Quanzhou 362000, P.R. China


<sup>b</sup> Department of Histology and Embryology, Hunan University of Medicine, Huaihua, P. R. China

<sup>c</sup> The Second Affiliated Hospital Of Fujian Medical University, Quanzhou 36200, P.R.China

<sup>d</sup> Xiamen Institute of Rare Earth Materials, Haixi Institute, Chinese Academy of Sciences, Xiamen 361021, China

E-mail: smilepearl@126.com (Z Zheng), aziliying@163.com (L Chen), songliang@fjirsm.ac.cn (L Song)

<sup>#</sup> These authors contributed equally to this work.

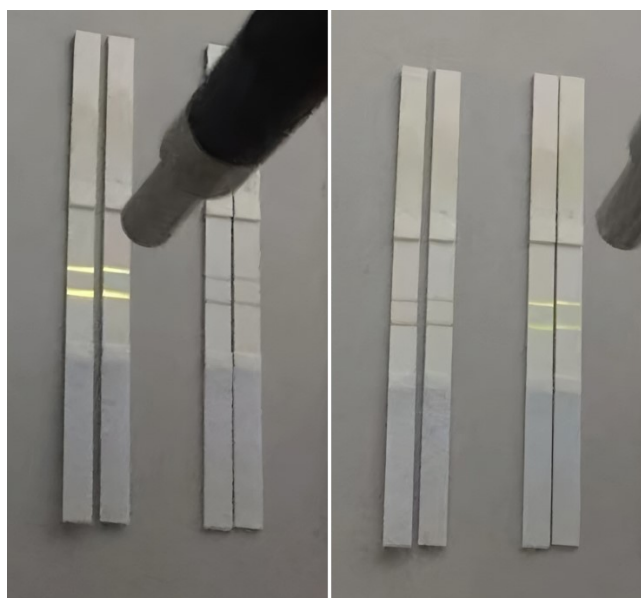

**Fig. S1** Photographs of upconversion fluorescence visualization detection results for FSH serum samples (concentration: 40.89 ng/mL (left) and 7.89 ng/mL (right)) under natural light condition.

**Table S1 Thermal stability test results of rare earth upconversion fluorescent immunochromatographic test strips**

	FSH		AMH	
Days	CV (%)			
	42.32 ng/mL	5.12 ng/mL	9.57 ng/mL	4.42 ng/mL
3	6.89	7.23	7.34	5.79
6	6.53	6.98	6.56	5.21
12	5.99	7.77	6.63	6.06
24	7.12	8.12	8.24	5.42
30	7.68	8.63	8.01	5.55